

## DESIGNED SPECIFICALLY FOR TEMPORARY HOSPITALS

The COV19 Response Bed is a manual medical bed that has been specifically designed to help manage the COVID 19 outbreak, and support the NHS in caring for and treating patients with the virus.

The COV19 Response Bed is completely manual so does not require power ensuring available power supply can be used for lifesaving medical equipment, such as, ventilators.

# COVID 19 RESPONSE BED




Designed specifically for use when temporary hospitals are setup to respond to an emergency situation that requires a functional bed.

The Field Response Bed is completely manual, it does not require a power supply, and ensures available power supply in Football Stadiums, exhibition halls, airport hangers etc can be used for lifesaving medical equipment.

## KEY FEATURES

- Manual Adjustable Backrest
  - o 30 degrees. Recommended during mechanical ventilation to reduce risk for ventilator-associated pneumonia.
  - o 30 to 45 degrees. Semi-Fowler position, used to promote oxygenation via maximum chest expansion and is implemented during events of respiratory distress
  - o 45 to 60 degrees. Fowler's position, facilitates the relaxing of tension of the abdominal muscles, allowing for improved breathing especially with respiratory complications.
  - o 70 degrees to facilitate feeding and improved postural drainage
- Low head end panel to facilitate intubation procedures
- Fixed height bed with folding legs for quick setup and storage
- No power supply requirement, suitable for makeshift areas
- IV Pole used to carry bags and deliver the necessary medications to patients
- Clamp on safety sides reduces the risk of falls

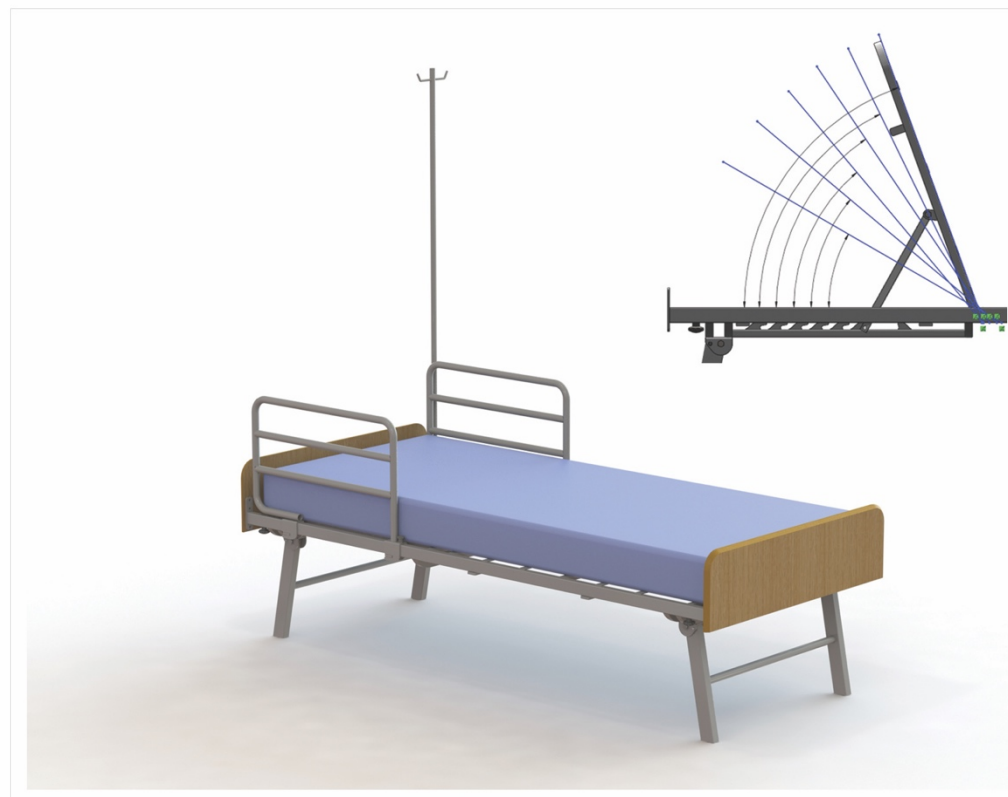

**Stock code**

COV19

**Description**

COVID 19 Response Bed

**Mattress platform dimensions**


2046mm long x 930mm wide

**Mattress platform height**

451mm

Enquire about this product - 01606 593 456 • [sales@renrayhealthcare.com](mailto:sales@renrayhealthcare.com)

### PRODUCT FEATURES

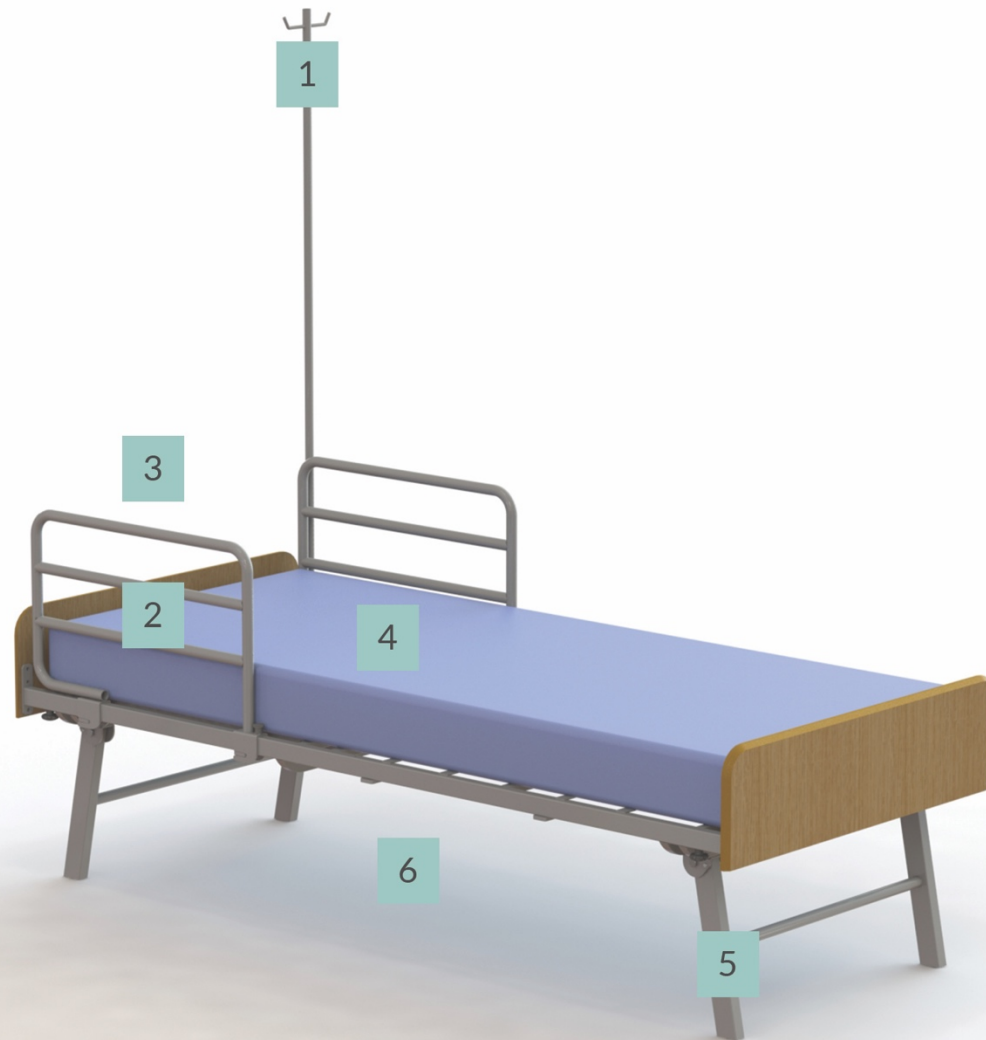

- 1 IV POLE** removeable, 1.5m fixed length
- 2 OR 4 PART SIDE RAILS** removeable, side rails can be fitted to both head and foot end as required
- 3 LOW HEAD END** for intubation access
- 4 BACKREST** multi-position backrest from 0 to 70 degrees for therapeutic positioning and feeding
- 5 CLEAR UNDER BED SPACE** provides easy access for cleaning and hoist access
- 6 STABLE FOLD AWAY LEGS** for compact storage. No castors, to prevent accidental bed movement in busy multi-bed ward areas or when patient exits the bed