


EN60601-2-52: Rail Ruler Guide

By Renray Healthcare

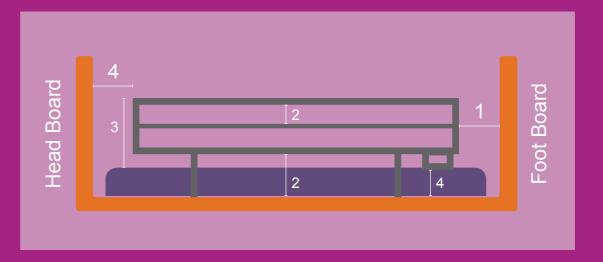

Renray Healthcare Limited, Road Five, Winsford Industrial Estate, Winsford, Cheshire CW7 3RB www.renrayhealthcare.com | E: sales@renrayhealthcare.com | T: 01606 593456 | F: 01606 861354

32cm

Measurement 1

Gap between foot board and end of side rail must be >318mm or <60mm

Crencay

Designed with you in mind y

As of the 1st April 2013 a new bed safety standard (EN60601-2-52:2009) came into force. This new Standard was introduced directly in response to a number of fatalities which had occurred with beds that meet the old bed standard (EN1970:2000). This new standard represents the current thinking in basic safety and essential performance of the electric profiling bed.

Measurement 3
Distance from top of mattress to top of side rails must be >220mm

22cm

12cm

Measurement 2
Gap between the side rail and mattress platform must be <120mm. Also applies to gaps in the side rails